

Table Of Contents

165-001 Blood Drawing Chair with Drawer Unpacking and Installation	1
Converting from Left-Hand Drawer Mounting to Right-Hand Drawer Mounting	1
Converting 165-002 Blood Drawing Chair with Fixed Arm from a Left-Hand Fixed Arm Mounting to a Right-Hand Fixed Arm Mounting	2
Height Adjustment for Rotating Arm Assembly, Fixed Arm Assembly, and Drawer Assembly	3
Rotating Arm Operation	3
Care and Maintenance	
Care of Upholstery	3
Care of Painted Metal Surfaces	3

MIDMARK®

Installation and Operation Instructions

Model 165-001 Blood Drawing Chair with Drawer
Model 165-002 Blood Drawing Chair with Fixed Arm

165-001 Blood Drawing Chair with Drawer Unpacking and Installation

1. Remove the **165-001 Blood Drawing Chair with Drawer** from the carton.
2. For packaging reasons, the drawer section of the chair is installed over the seat section of the chair. Loosen knob "1" (Fig. 1) and remove the drawer section from the "H"-frame (2, Fig. 1). Rotate the drawer section 180°. Insert the drawer section back into the "H"-frame and tighten knob "1" (Fig. 1).

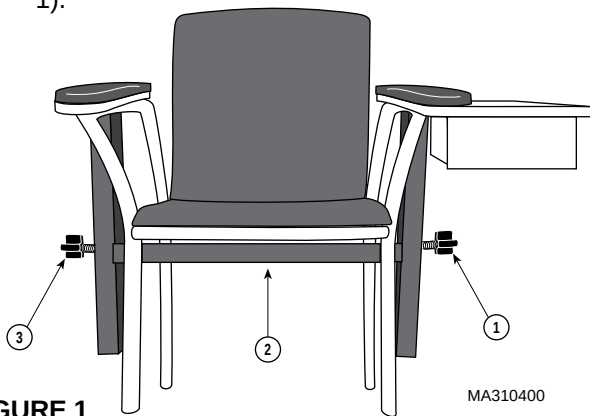


FIGURE 1

MA310400


Converting from Left-Hand Drawer Mounting to Right-Hand Drawer Mounting

NOTE

All references to left and right are from the patient's position while seated in the chair.

1. The **165-001 Blood Drawing Chair with Drawer** is shipped from the factory with the drawer assembly mounted on the left-hand side. To change the drawer assembly to a right-hand mounting, loosen knobs "1" and "3" (Fig. 1).
2. Remove the drawer from the "H"-Frame and set it aside. Remove the rotating arm from the "H"-Frame and install it into the other side of the "H"-Frame. Be sure that the locking-knob (2, Fig. 2) of the rotating arm (1, Fig. 2) is positioned to the front of the chair.

3. Take the drawer assembly (removed in Step 2) and place the white Formica top side down. Be sure to place the white Formica top on a soft surface so that it won't be scratched.

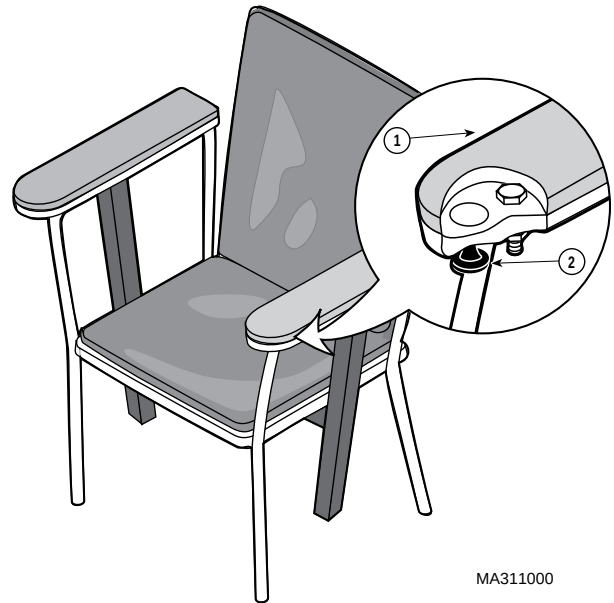

MA311000

FIGURE 2

4. Using a Phillips screwdriver, remove the four wood screws (1, Fig. 3) from the Formica top.
5. Rotate the drawer and the drawer slide 180°. Align the four holes in the Formica top with the four holes in the drawer slides. Reinstall the four Phillips head wood screws that were removed.

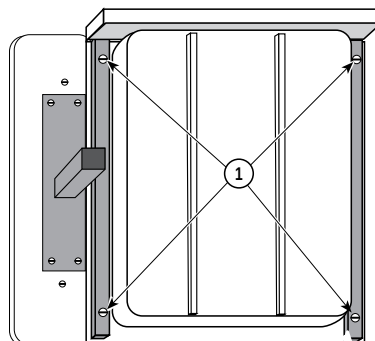


FIGURE 3

MA310502

6. Tighten the four screws with a Phillips screwdriver.

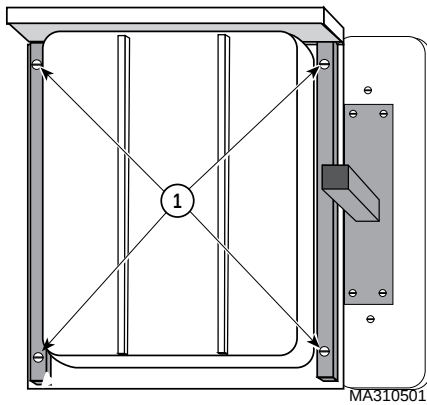


FIGURE 4

7. Place the drawer assembly into the remaining side of the "H"-Frame. Tighten knobs "1" and "2" (Fig. 5).

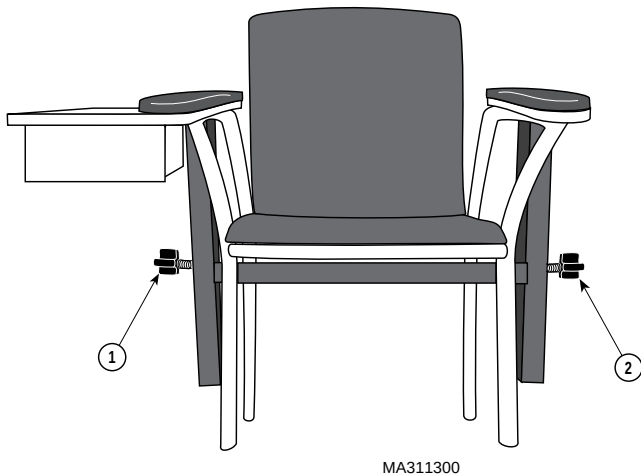


FIGURE 5

Converting 165-002 Blood Drawing Chair with Fixed Arm from a Left-Hand Fixed Arm Mounting to a Right-Hand Fixed Arm Mounting

1. The **165-002 Blood Drawing Chair with Fixed Arm** is shipped from the factory with the fixed arm assembly (1, Fig. 6) mounted on the left-hand side. To convert the fixed arm assembly to a right-hand mounting, loosen knobs "2" and "3" (Fig. 6).

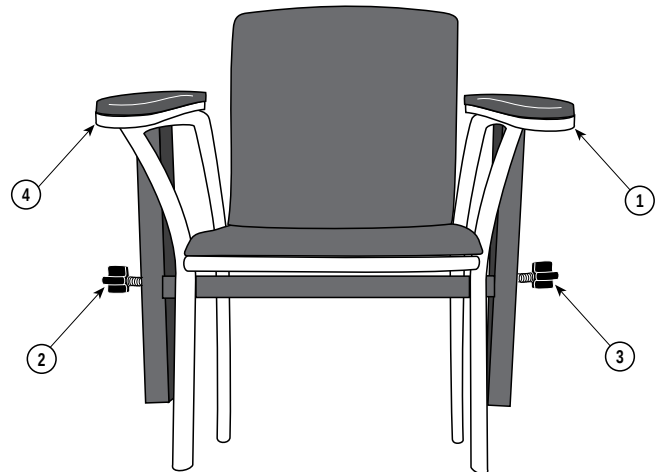


FIGURE 6

2. Remove the fixed arm assembly (1, Fig. 6) on the left-hand side from the "H"-Frame and set the arm assembly aside.
3. Remove the rotating arm assembly (4, Fig. 6) from the "H"-Frame and install it into the other side of the "H"-Frame. Be sure the locking knob (2, Fig. 2) of the rotating arm (1, Fig. 2) is positioned to the front of the chair.
4. Take the fixed arm assembly (removed in Step 2) and install it into the remaining side of the "H"-Frame. Tighten knobs "1" and "2" (Fig. 7).

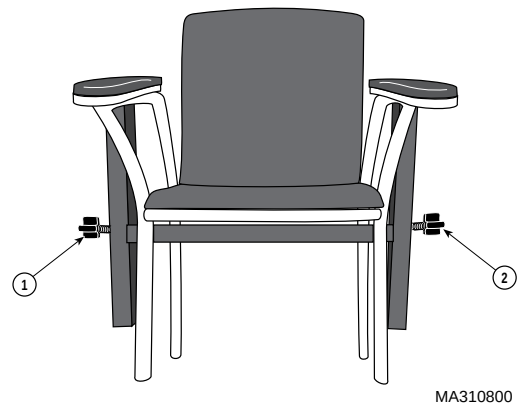


FIGURE 7

Height Adjustment for Rotating Arm Assembly, Fixed Arm Assembly, and Drawer Assembly

The height of the rotating arm assembly, fixed arm assembly, and the drawer assembly can be adjusted up or down for patient comfort. This can be accomplished by the following steps:

1. Loosen the knob (1, Fig. 8) on the "H"-Frame (2, Fig. 8).
2. Raise or lower the assembly as required.
3. Tighten the knob.

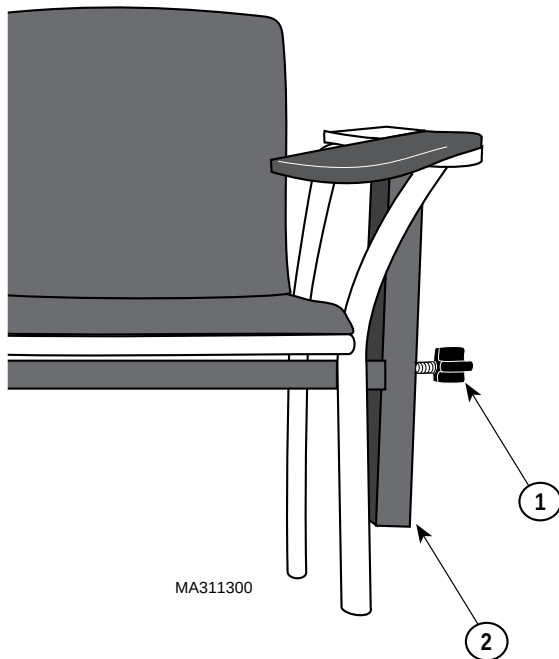


FIGURE 8

Rotating Arm Operation

The pad on the rotating arm assembly can be rotated and locked at any angular position. This can be accomplished by the following steps.

1. Loosen knob "2" (Fig. 2) below the arm pad.
2. Position the pad as required.
3. Tighten knob "2" (Fig 2).

Care and Maintenance

Care of Upholstery

The upholstery material used to cover the mattress is resistant to most medicinal-type stains, but may be damaged by solvents and dyes. Regular care should be maintained by daily wiping with a damp cloth or sponge and periodic cleaning with a mild soap and water solution.

Any fluid spilled on the upholstery should be removed as quickly as possible.

In the case of stain, attempt to remove it with soap and water.

Care of Painted Metal Surfaces

All painted metal surfaces should be wiped with a clean soft cloth weekly and periodic applications of paste wax made to all surfaces to preserve the finished luster.



Midmark Corporation, Versailles, Ohio 45380 U.S.A.
937-526-3662 TELEX 288007 FAX 937-526-5542